



Surgical Apparels

A healthcare professional, likely a nurse or surgeon, is shown from the chest up, wearing a blue surgical cap, a light blue surgical mask, and white gloves. They are looking upwards and to the right. The background is a blurred clinical setting with yellow and white curved shapes overlaid on the image.

ABOUT US

Welcome to 360° Infection Prevention Solutions Where Protection Comes Full Circle.

In the world of healthcare, safety is not a choice — it's a commitment. At 360° Infection Prevention Solutions, we've made that commitment to our mission by delivering a comprehensive suite of solutions that protect every touchpoint of surgical and clinical care. From the precision of our **Surgical Apparels** to the power of our **Disinfection & Sterilization** range, from expert-crafted **Hand Hygiene and Skin Preparation** products to the latest in **General and Advanced Wound Care** - we leave no corner uncovered.

We don't stop there. Our specialized offerings extend to **Ostomy Care** solutions that restore dignity and comfort, and **Sterilizers with Accessories** that ensure every tool is ready for the next life-saving procedure.

At 360°, we don't just provide products — we create a circle of trust, safety, and innovation around patients, practitioners, and healthcare spaces. Every solution we deliver is thoughtfully designed to drive better hygiene, stronger safety standards, and superior patient outcomes.

Because when it comes to infection prevention, every degree matters.

TABLE OF CONTENT

01.	Sub Vertical Overview	-----	3
02.	Surgical Apparels	-----	4-10
03.	Surgical Caps	-----	11-12
04.	Surgical Mask	-----	13
05.	Surgical Drapes	-----	14-21
06.	Upcoming Drapes	-----	22
07.	Mayya Kits	-----	22



ZERO COMPROMISE

01

Surgical Apparels

SURGICAL APPAREL

360° Infection Prevention Solutions | Surgical Apparel

At Meril, our **Infection Prevention & Control Division** is dedicated to advancing infection prevention in surgical settings through our specialized **360° Infection Prevention Solutions**. Our Surgical Apparel is purpose-built to raise the standards of protection in **Operating Theatres (OTs)** and other critical care areas.

We are committed to mitigating infection risks associated with surgeries by delivering a robust, clinically validated portfolio that empowers healthcare professionals with **protection, confidence, and compliance** with **global safety standards**.

OUR PRODUCT PORTFOLIO

➤ MAIRA™ Surgical Gowns

Superior barrier protection with certified **AAMI Level 1 to 4** standards, ensuring optimal **safety, comfort, and breathability** for extended surgical procedures.

➤ MEVEL™ Surgical Drapes


Precision-engineered sterile draping systems are designed to provide **comprehensive coverage and reliability** across a wide range of surgical interventions.

➤ MAYYA™ Surgery-Specific Kits

Procedure-customized surgical kits that streamline workflows, enhance surgical efficiency, and reduce the risk of contamination through standardized, high-quality components.

➤ MAIRA™ Surgical Caps

➤ MAIRA™ Surgical Mask

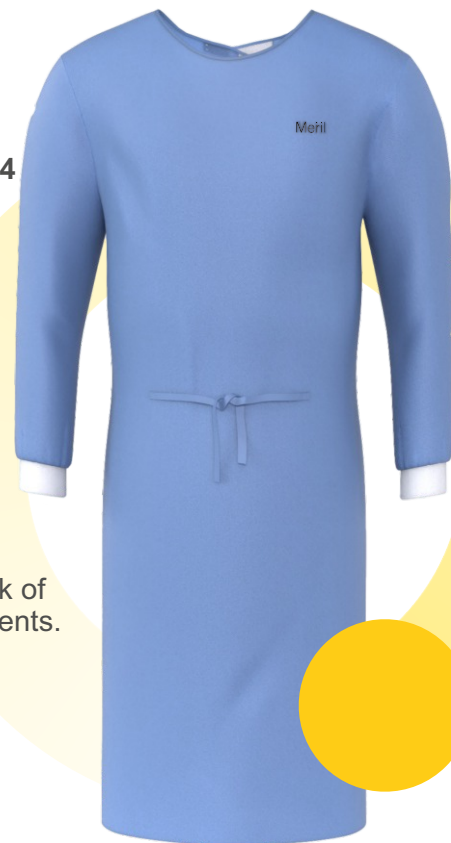

MAIRA™ SURGICAL GOWN


**Advanced Protection. Trusted Comfort.
Surgical Confidence.**



**High Performance
Breathable Gown**



**Meets AAMI
level Standards**

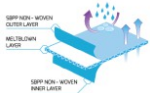

**Made of Polypropylene
SSMMS Fabric Hydrophobic**



ULTRASONIC TIE



KNITTED CUFF



**ULTRASONIC
SEAM**



**COLLAR VELCRO
DESIGN**

ENHANCED BARRIER PROTECTION

- **Breathable fabric** that resists **fluid penetration and blood strike-through**
- Effective **microbial penetration resistance** for infection Prevention
- **Ultrasonically sealed seams** for enhanced durability and fluid protection

SAFE FOR SKIN

- Certified **free from skin irritation, sensitization, and cytotoxicity** for patient and user safety

CUSTOM FIT

- Equipped with **Velcro hook and loop** closures for an adjustable neckline to ensure a personalized fit.

VIRTUALLY LINT-FREE


➤ Breathability

Advanced moisture vapor transmission allows for breathability compared to standard SMS gowns, keeping clinicians cool and comfortable

Sterility & Packaging

- **EO Sterilized** as per **ISO 11135 standards**
- **Book folded** for **aseptic handling and easy usage**
- **Pouch is packed** in **certified and validated medical-grade packaging**

FEATURES

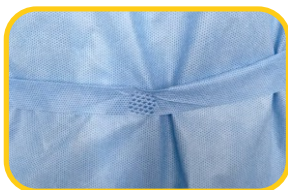

ULTRASONIC TIE



ULTRASONIC SEAM



WATER REPELLANT



KNITTED CUFF



COLLAR VELCRO DESIGN



MEETS AAMI LEVEL STANDARD

MAIRA™ Surgical Gown Variants

MAIRA™ ISOLATION GOWN

► Lightweight Protection. Everyday Confidence.

MAIRA™

Isolation Gowns are designed to provide essential protection for healthcare professionals in low-risk environments.

► Standard Specification

Type – Isolation Gown

Material - 35 GSM SMS Medical Textiles

Color - Blue

AAMI Level - 1

Design - Straight Shoulders - Velcro Neck

Units Per Case – 10 Pcs

Sterilization – EO Sterile



Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

MAIRA™ STANDARD GOWN

► Reliable Protection. Everyday Performance.

MAIRA™

Standard Gowns are designed to meet the needs of healthcare professionals in procedures involving low to moderate fluid exposure.

► Standard Specification

Type – Standard Gown

Material – 43 GSM, SMS Medical Textiles

Color – Blue

AAMI Level – 2

Design - Straight Shoulders - Velcro Neck

Cuff – Elastic Cuff

Sterilization – EO Sterile

Units Per Case – 10 Pcs



Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

MAIRA™ REINFORCED GOWNS

➤ Extra Protection Where It Matters Most.

MAIRA™

Reinforced Gowns are designed for procedures involving moderate fluid exposure, offering AAMI Level 3 protection. The gowns combine strength, comfort, and targeted safety to meet the needs of high-demand clinical settings.



➤ Standard Specification

Type - Reinforced Gowns

Material - 43 GSM, SMS Medical-grade textiles,
Reinforced in Critical Zone

Color - Blue

AAMI Level - 3

Design - Straight Shoulders - Velcro Neck

Cuff - Elastic Cuff

Sterilization - EO Sterile

Units Per Case - 10 Pcs

Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

MAIRA™ FULLY REINFORCED GOWN

➤ Defining Protection. Delivering Assurance.

MAIRA™

Fully Reinforced Gowns set a new benchmark in surgical apparel protection. Built on a 43 GSM SMS base fabric, these gowns are entirely reinforced from collar to hemline and across full sleeves with a **30 GSM highly breathable film**, delivering AAMI Level 4 protection.



➤ Standard Specification

Type - Fully Reinforced Gowns

Material - 43 GSM, SMS Medical-grade textiles,
Reinforced in Critical Zone (**Sleeves, Chest, and Arm Pit**)

Color - Blue

AAMI Level - 4

Design - Straight Shoulders - Velcro Neck

Cuff - Elastic Cuff

Sterilization - EO Sterile

Ideal Applications - Cardiovascular, Orthopaedic,
high fluid exposure Cases

Units Per Case - 10 Pcs

Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

MAIRA™ PREMIUM SURGICAL GOWN

► Ultimate Protection. Professional Confidence.

MAIRA™

The MAIRA™ Premium is built for the highest-risk surgical environments. With 60 GSM trilaminate fabric and AAMI Level 4 certification, this gowns provide maximum fluid and pathogen protection without compromising comfort.

► Standard Specification

Type - Pro Surgical Premiere Gowns

Material - 60 GSM, trilaminate Medical-grade textiles

Color - Blue

AAMI Level - 4

Design - Straight Shoulders - Velcro Neck

Cuff - Elastic Cuff

Sterilization - EO Sterile

Ideal Applications - Cardiovascular, Orthopaedic, high fluid exposure Cases

Units Per Case - 10 Pcs



Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

MAIRA™ SUPREME GOWN

► The Right Balance Between Comfort & Safety.

MAIRA™

Supreme Gowns are designed with 51 GSM advanced fabric, delivering the highest protection (AAMI Level 4) while ensuring all-day comfort. Built with **Breathable Viral Barrier (BVB)** film technology.

► Standard Specification

Type - Breathable / BVB Supreme Gown

Material - 51 GSM, SMS Medical-grade textiles, with Breathable Viral Barrier (BVB) film technology

Color - Blue

AAMI Level - 4

Design - Straight Shoulders - Velcro Neck

Cuff - Elastic Cuff

Sterilization - EO Sterile

Ideal Applications - Cardiovascular, Orthopaedic, high fluid exposure Cases

Unit Per Case - 10 Pcs




Size	Length	Width
L	124 cm	160 cm
XL	134 cm	160 cm
XXL	144 cm	160 cm
XXXL	150 cm	160 cm

UNDERSTAND THE AAMI LEVEL GUIDELINES

ANSI/AAMI PB70:2012 -Liquid Barrier Performance.

The AAMI Level Standards define the barrier performance requirements of surgical gowns, drapes, and other protective apparel. These standards enable healthcare providers to select the appropriate level of protection based on the expected fluid exposure and associated infection risk of medical procedures.

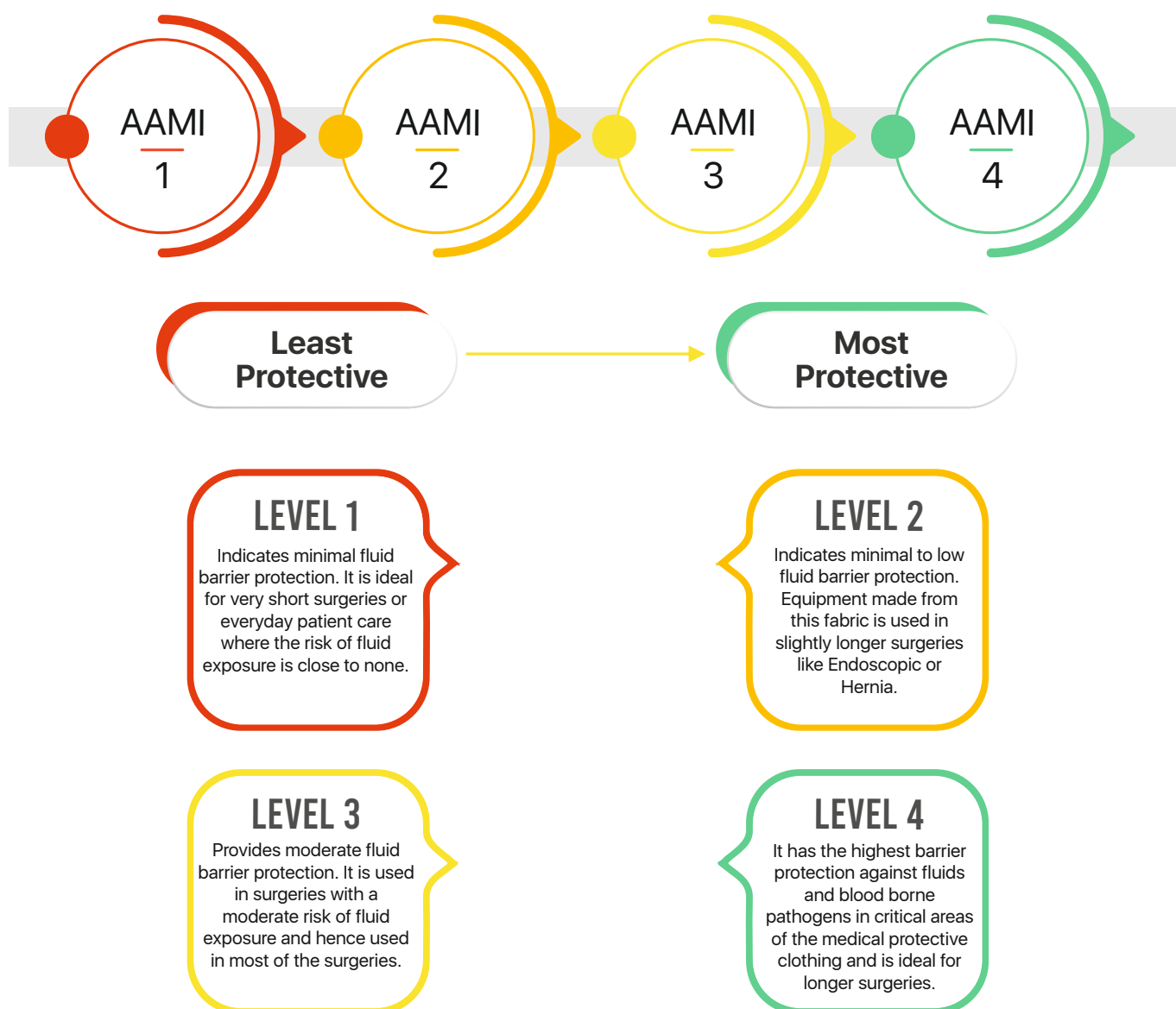

{***AAMI** - ASSOCIATION FOR THE ADVANCEMENT OF MEDICAL INSTRUMENTATION}

{***ANSI** - AMERICAN NATIONAL STANDARDS INSTITUTE}

AAMI Level	Protection Level	Typical Use	Fluid Resistance	Test performed
Level 1	Minimal Protection	For general patient care or simple examination procedures, where minimal fluid exposure is anticipated	Limited resistance to fluids like water or blood	Impact Penetration Test (ISO 22610) Barrier Protection Test (ASTM F1670)
Level 2	Low to Moderate Protection	Suitable for procedures where there is a low to moderate risk of fluid exposure, such as during blood draws, suturing, or minor surgeries	Moderate resistance to fluids and penetration	Impact Penetration Test (ISO 22610) Barrier Protection Test (ASTM F1670)
Level 3	Moderate to High Protection	Used for moderate to high-risk surgeries or procedures, such as orthopedic surgeries, catheterization, and diagnostic testing, where moderate fluid exposure is expected	High resistance to fluid splashes and penetration	Impact Penetration Test (ISO 22610) Barrier Protection Test (ASTM F1670) Water Pressure Test (AAMI PB70) Blood Penetration Test (ASTM F1670)
Level 4	High Protection	For high-risk procedures involving high fluid exposure, such as trauma surgeries, organ transplants, and dealing with infectious disease patients, where the risk of contamination is severe	Maximum resistance to fluid penetration, ensuring complete protection in critical settings	Impact Penetration Test (ISO 22610) Barrier Protection Test (ASTM F1670) Water Pressure Test (AAMI PB70) Blood Penetration Test (ASTM F1670) Cyclic Testing (AAMI/ANSI Pb70)



MAIRA™ SURGICAL CAP

Secure. Sterile. Comfortable.



MAIRA™ Surgical Cap is a high-quality head covering designed for use in sterile surgical environments. It ensures comprehensive protection by containing hair and reducing microbial dispersion from the scalp and hair, thus contributing to a sterile field.

Feature	Description
Material	SMS (Spunbond-Meltblown-Spunbond) Non-woven Fabric
Weight	25-30 GSM
Design	Elasticated band for snug fit
Breathability	Air-permeable, minimizes sweating
Barrier Protection	Limits hair shedding and microbial contamination
Comfort Fit	Soft texture, skin-friendly
One Size Fits All	Universal fit with secure hold
Latex-Free	Reduces allergy risk
Colour Options	Standard Blue, White (custom colours optional)

BENEFITS

- **Infection Prevention:** Prevents hair and scalp microbes from contaminating the sterile field.
- **Enhanced Comfort**
- **Secure Fit**
- **Hypoallergenic**



MAIRA™
Surgical Cap

Pack Size

Available in a Pack of 50 Pcs

Sterilization – EO Sterilized

MAIRA™ BOUFFANT CAP

Full Coverage, Maximum Protection, Ultimate Comfort.



The **MAIRA™ Bouffant Cap** is a soft, lightweight, and breathable head covering designed to fully encapsulate hair, reducing the risk of hair and microbial shedding in sterile and clean environments. Its roomy design accommodates various hair volumes, ensuring both comfort and infection Prevention.

Feature	Description
Material	Spunbond Polypropylene Non-woven Fabric
GSM Range	10-15 GSM (customizable)
Design	Pleated, Full-Coverage Circular Design
Elastic Band	Soft, non-binding elastic for secure yet gentle fit
Air Permeability	High breathability to reduce sweating
Universal Fit	Available in 18", 21", and 24" sizes
Latex-Free	Prevents latex allergy concerns
Color Options	Blue

BENEFITS

- **Enhanced Infection Prevention:** Prevents hair and particulate contamination in sterile zones
- **Comfort for Long Wear**
- **Versatile Sizing**
- **Safe for All Users**
- **Economical**



MAIRA™
Bouffant Cap

Pack Size

Pack of 50 Bouffant Caps

MAIRA™ SURGICAL MASK

Breathe Easy. Stay Protected.



The **MAIRA™ Surgical Mask** is carefully engineered to provide effective respiratory protection for healthcare personnel, patients, and the general public. Designed to cover the nose and mouth, it filters airborne pathogens while allowing for comfortable, normal breathing. Suitable for use in the Operating Theatre (OT), isolation settings, and public health emergencies.

Feature	Description
Protection	Protects against pathogenic biological airborne particulates
Layer Construction	3-Layer Protection: Outer Hydrophobic Layer, Middle Meltblown Filter Layer, Inner Soft Absorbent Layer
Filtration Efficiency	≥98% Bacterial Filtration Efficiency (BFE)
Breathability	Low breathing resistance for comfort
Design	Pleated design for face coverage and easy adjustment
Attachment Type	Strong, soft tie attachments for a secure and adjustable fit
Nose Clip	Integrated, flexible nose piece for improved seal
Comfort	Soft inner lining is gentle on the skin
Usage	Disposable, single-use
Standard Compliance	Follows ASTM Level 2 / EN14683 Type IIR equivalent standards

BENEFITS

- **Optimal Protection:** Filters airborne bacteria and particulates with high efficiency
- **Comfortable Fit**
- **Skin-Friendly**
- **Enhanced Breathability**
- **Versatile Applications**



MAIRA™
Surgical Mask

Pack Size

Pack of 50

MEVEL™ SURGICAL DRAPE

Precision. Protection. Performance.


Key Benefits of Our Drapes

- 100% Impermeable — Complete protection against fluids
- High Absorption Capacity & Speed — Rapid absorption minimizes pooling
- Low Linting — Maintains a sterile field by reducing particle contamination
- Soft Feel & Body Conformity — Ensures patient comfort and secure draping
- Latex-Free — Safe for latex-sensitive users and patients

ENHANCED FEATURES

- **Royal Blue SSMMS Non Woven Drape Fabric:**
Reduces glare from operating lights, improving visibility
- **Standardized Reinforcement:**
Critical zones reinforced with advanced materials for added protection
- **Clear Visual Contrast:**
Differentiates drape fabric from reinforced areas for better orientation during procedures
- **Enhanced Sizing Options:**
Expanded reinforcement sizing on select codes for optimal fit and coverage

{*SSMMS - Spunbond Spunbond Meltblown Meltblown Spunbond}

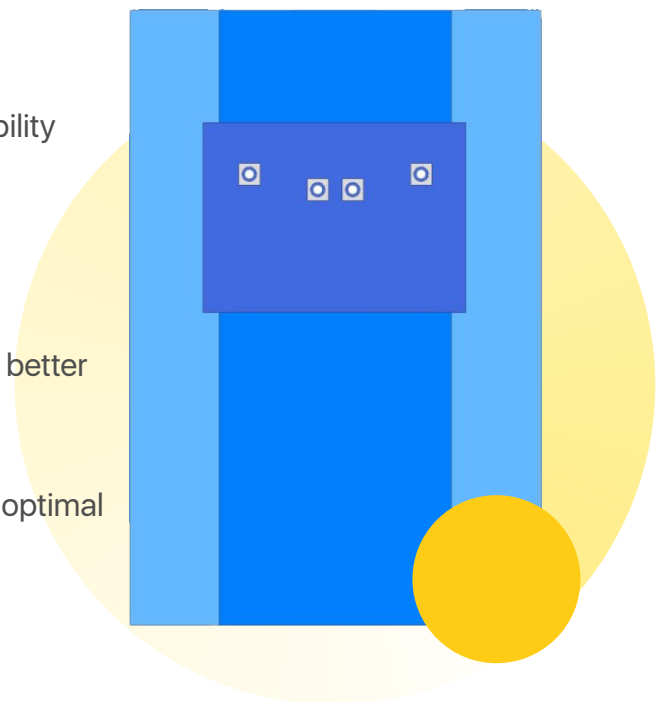

Mevel

MEVEL™ ORTHOPEDIC DRAPE


Mevel™ Orthopedic Drapes are made from SSMMS non-woven fabric combined with absorbent non-woven layers. These drapes are designed to effectively isolate the surgical site from other areas of the patient's body and the non-sterile zones of the OR table, thereby helping to reduce the risk of surgical site infections (SSI). **Resistant to fluid penetration**, fluid concentration, and blood strike-through, **Free from skin irritation, sensitization, and cytotoxicity**, Resistant to **microbial penetration**, **Virtually lint-free** and **abrasion-resistant**. Drapes are **pouch-packed in certified and validated packs**, **Ethylene Oxide (EO) sterilized** as per **ISO 11135** standards.

1 MEVEL™ HIP U DRAPE

The **Mevel Hip U Drape** is a high-performance surgical drape designed for a wide range of procedures that demand efficient fluid management and optimal site access.

Features :-

- Centrally located **U-shaped split measuring 10 cm x 80 cm**
- Drape is reinforced with premium SMPE absorbent fabric
- Equipped with **Tube Holder** and **Cord Organizer**

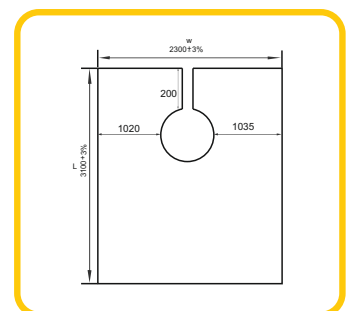


Size: 310 cm × 230 cm

2) MEVEL™ KNEE O DRAPE

The **Mevel Knee O Drape**, measuring 300 × 230 cm, is a specialized surgical drape designed for knee and lower limb procedures requiring wide coverage and enhanced fluid control.

Features :-

- **Central hole with a 2.5-inch diameter**
- **Layer of latex-free material**
- Drape is **reinforced with premium SMPE absorbent fabric** to ensure high fluid absorbency
- Equipped with **Tube Holder** and **Cord Organizer**

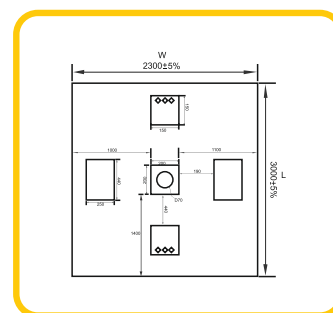


Size: 300 cm × 230 cm

3) MEVEL™ EXTREMITY DRAPE

The Mevel Extremity Drape, measuring 360 cm × 240 cm, is designed to provide comprehensive coverage and effective fluid management for upper and lower extremity surgical procedures.

Features :-

- Central **2.5-inch** diameter hole that allows precise access to the surgical site
- **Latex free material**
- Drape is **reinforced with premium SMPE absorbent fabric** to ensure high fluid absorbency
- Equipped with **Tube Holder** and **Cord Organizer**

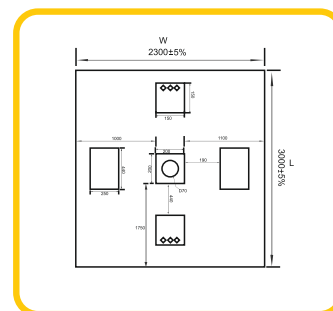

Size: 360 cm × 240 cm





MEVEL™ CARDIOLOGY DRAPE

Mevel Angiography drapes, made from SSMMS non-woven fabric and absorbent non-woven fabric, help isolate the surgical site from other areas of the patient's body and non-sterile areas of the OR table, contributing to a reduction in the risk of surgical site infection (SSI). The fabrics meet AAMI PB 70 level standards.


The fabric used in making the drapes are resistant to fluid penetration, fluid concentration, and blood strike-through; they are free from skin irritation, skin sensitization, and cytotoxicity. Also resist microbial penetration, are virtually lint-free, and are abrasion-resistant.

1 MEVEL™ FEMORAL ANGIOGRAPHY DRAPE

Mevel Femoral Angiography Drape Ultra of size 300 cm × 230 cm,

Features :-

- Two femoral fenestration of 11 cm diameter
- Drape is **reinforced with premium SMPE absorbent fabric** to ensure high fluid absorbency

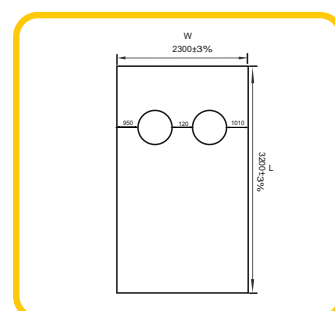

Size: 300 cm × 230 cm

MEVEL™ ENT / EYE / DENTAL DRAPE


Mevel

1 MEVEL™ DENTAL DRAPE

The **Mevel Dental Drape**, measuring **50 cm × 80 cm**, is specifically designed for dental and maxillofacial procedures

Features :-

- Centrally placed **10 cm × 10 cm fenestration** with pouch
- Drape is **reinforced with premium SMPE absorbent fabric** to ensure high fluid absorbency

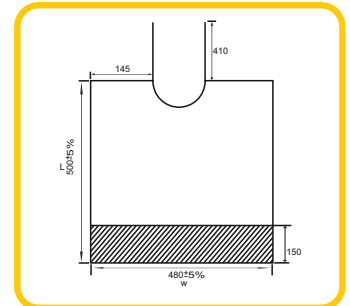


Size: 50 cm × 80 cm

2 MEVEL™ EYE DRAPE

Mevel Eye Drape, sized **80 cm × 90 cm**, is designed to provide optimal coverage and fluid control during eye surgeries and ophthalmic procedures.

Features :-

- **8 cm × 8 cm fenestration** with pouch
- Reinforced with an **SMPE Absorbent** layer to manage fluids effectively
- Transparent incise film that ensures secure adhesion to patient skin

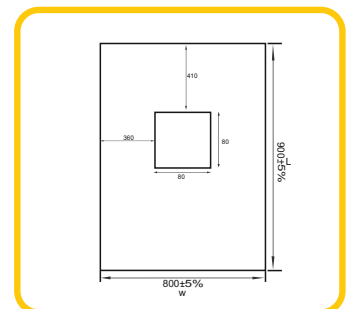

Size: 80 cm × 90 cm

{*SMPE Spunbond-Meltblown-Polyethylene-Laminated





MEVEL™ GYNECOLOGY DRAPES


Mevel Gynecology Drapes are made from SSMMS non-woven fabric combined with absorbent non-woven layers. These drapes are designed to effectively isolate the surgical site from other areas of the patient's body and the non-sterile zones of the OR table, thereby helping to reduce the risk of surgical site infections (SSI). **Resistant to fluid penetration**, fluid concentration, and blood strike-through, **Free from skin irritation, sensitization, and cytotoxicity**, Resistant to **microbial penetration**, **Virtually lint-free** and **abrasion-resistant**. Drapes are **pouch-packed in certified and validated packs**, **Ethylene Oxide (EO) sterilized** as per **ISO 11135** standards.

1 MEVEL™ C SECTION DRAPE

The Mevel C-Section Drape with Fluid Collection Pouch, measuring 190 cm × 230 cm, is specially designed for cesarean section procedures, offering comprehensive coverage, fluid control, and safety.

Features :-


- Large 20 (diameter) cm fenestration
- Reinforced backdrop to provide enhanced protection against fluid penetration and pooling
- Integrated fluid collection pouch
- Equipped with **Tube Holder** and **Cord Organizer**

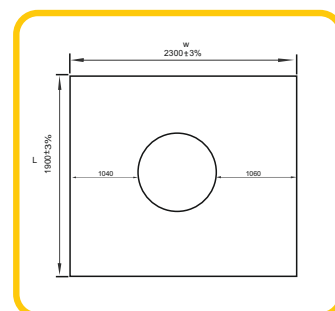

Size: 190 cm × 230 cm


MEVEL™

SURGICAL DRAPE SHEET

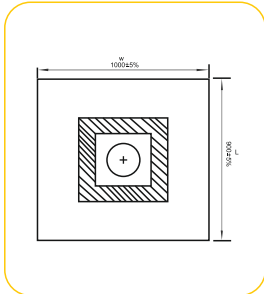


Mevel™ Fenestrated Drape

- Diameter of fenestrations – 11 Cm
- Size : 90 cm - 100 cm

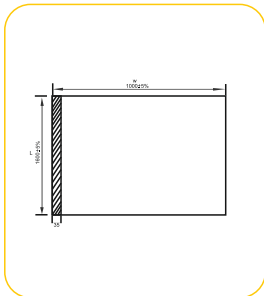


Mevel™ Medium Drape

- With adhesive absorbent
- Size : 160 cm - 100 cm

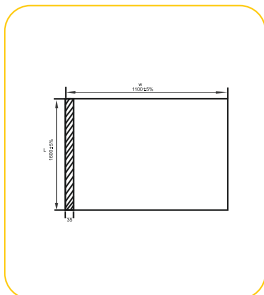


Mevel™ Half Drape

- With Adhesive
- Size : 160 cm - 110 cm

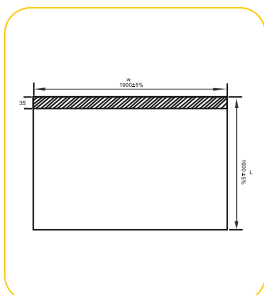


Mevel™ Large Drape

- With adhesive
- Size : 190 cm - 160 cm

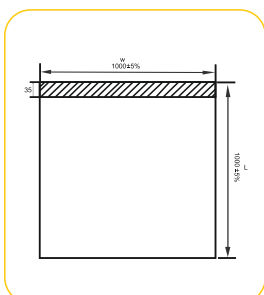


Mevel™ Utility Drape

- With adhesive
- Size : 100 cm - 100 cm

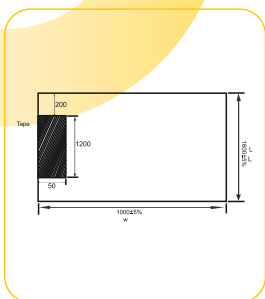


Mevel™ Top Drape

- With Adhesive
- Size : 100 cm - 160 cm

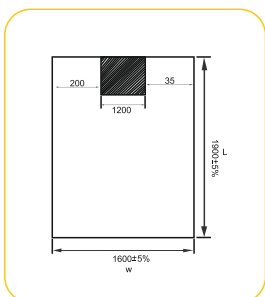

Mevel™ Bottom Drape

- with Adhesive
- Size : 190 cm - 160 cm

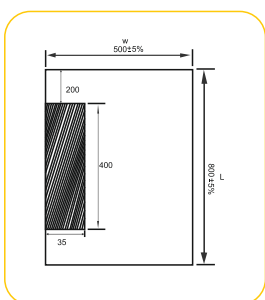

Mevel™ Side Drape

- With Adhesive
- Size : 50 cm - 80 cm



UPCOMING DRAPES

Category	Type
ORTHOPAEDIC DRAPES	Bilateral Drape Arthroscopy Drape Lamino Spinal Drape Endoplasty Spine Drape
Category	Type
NEUROLOGY DRAPES	Craniotomy Drape
Category	Type
UROLOGY DRAPES	TURP Drape PCNL Drape
Category	Type
CARDIOVASCULAR DRAPES	CABG Drape
Category	Type
GENERAL SURGERY DRAPES	Laparoscopy Drapes Laparotomy Drapes Adhesive Drapes

THE STANDARDIZED KITS AND CUSTOMIZED SURGERY-SPECIFIC KITS WILL BE AVAILABLE SOON.

THE NINE STEPS OF HOW TO PUT ON A MAIRA SURGICAL GOWN ASEPTICALLY

1.

Grasp the sterile gown and step into an area where the gown may be opened without risk of contamination. Put both hands into armholes. Hold the gown away from the body, and allow it to unfold completely.



2.

Slip both hands and forearms into armholes and sleeves; keep hands at shoulder level and away from the body.



3.

If you are open glove donning, pull cuff to thumb level. If you are closed glove donning, push the hands to the edge of the cuff.



4.

The circulating nurse should pull the gown over shoulders, touching only the inside of the gown.



5.

The circulating nurse hook and loop fixation.





should arrange



Retrieve the outer belt by pulling it out of the belt card held by the circulating nurse and tie it together with the other outer belt on the left-hand side.

9.

8.

Hand the belt card to the circulating nurse, then make a three-fourth turn to the left while the circulating nurse extends the outer belt to its full length.



Then don gloves. After gloving, grasp the belt card with both hands, separate the left outer belt from the belt card and keep the belt in the left hand.

7.

6.

The circulating nurse should tie the inner belt of the gown.



35 Meril

Subsidiaries

12 Academies

Globally

10000+
EMPLOYEES

250+
PRODUCTS

158+
COUNTRIES



 **Manufactured by:**

Meril Medical Innovation Pvt. Ltd.

Ground & First Floor, Survey No. 173/4 &
First Floor, H1-H3, Meril Park,
Survey No. 135/2/B & 174/2, Muktanand Marg,
Chala, Vapi - 396 191, Gujarat, India

Customer Care Contact:

Email: askinfo@merillife.com

Marketed By:

Meril Life Sciences India Pvt. Ltd.

6th Floor, B wing, Bonanza Building,
Sahar Plaza, Chakala JB Nagar,
Nr Kohinoor Hotel, Nr Chakala Metro Station,
Andheri East. 400059